




JAGA penny Receiver options

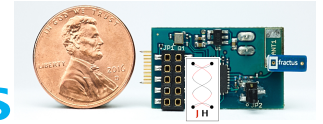

2018

JAGA penny is a very low power wireless neural recording solution. The signal range is around 5m.

We provide two versions of receivers (A, B). The option B comes with the USB cable (See Right). (The option C phased out.) Users can test these two options and choose the receiver that best suits their experiments.

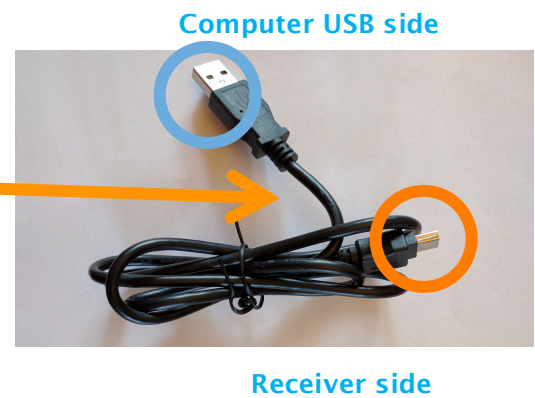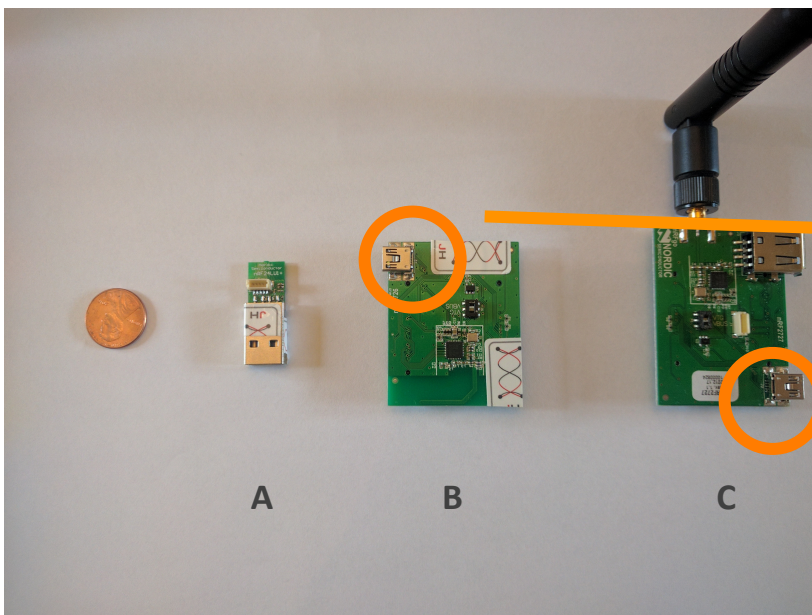

Receiver options:

- **A:** Most elegant and pretty option. It can directly be plugged into your computer. With this option, your computer needs to be very close to the receiver. The range can be limited.
- **B:** The receiver and the computer can be further away with the receiver and the computer connected to the USB cable shown on the right hand side. The plug for the receiver side is shown with a yellow circle, and the plug for the computer size is shown with a blue circle.
- **C:** With an external antenna, it can boost up the signal strength. We phased out of this option. We found that the performance was not any better than A, B.
- **Please plug the receiver first into the computer *before* power the JAGA Penny device.**